

Submittal Data Sheet

Project Information

Project Number _____ Approval _____



64-0603



64-0604-6B




64-0003 with 64-0605




64-0603-6B with 64-0606 mounted on slide bracket with suction regulator



64-0606A

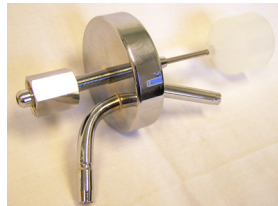

64-06CAP

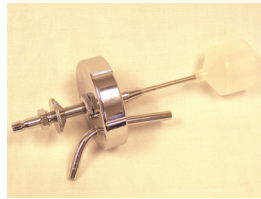

64-06CAP-6B



64-0603 with 64-0605 mounted on slide bracket with suction regulator

Tri-Tech Medical fluid waste collection bottles are designed with safety and economy in mind. Polycarbonate bottles are durable and reusable. Polycarbonate bottles may be autoclaved or sterilized with sterilizing gases and re-used many times. Polycarbonate bottles permit easy visual inspection of the collection bottle contents.

Collection bottle and cap assemblies are designed for high air flow capability. The assemblies consist of the polycarbonate bottle and metal cap assembly. The cap assemblies are designed with float shut-off, vacuum tube extension, and cap. The cap assemblies include a DISS 1220 female wrench tight nut and nipple top connection or a suction tubing nipple connection. The cap assemblies with suction tubing nipple connections feature locking stems and nuts for use with suction bottle and cap assembly mounting brackets. The mounting brackets are designed to mount in slide brackets.

Part Number	Description
64-0003	1200 cc (1/2 gallon) Polycarbonate bottle only – packaged in cases of 6 bottles each
64-0004	2800 cc (1 gallon) Polycarbonate bottle only – packaged in cases of 6 bottles each
64-0603-6B	1200 cc (1/2 gallon) Polycarbonate bottle and cap assembly with suction tubing connection
64-0604-6B	2800 cc (1 gallon) Polycarbonate bottle and cap assembly with suction tubing connection
64-0603	1200 cc (1/2 gallon) Polycarbonate bottle and cap assembly with DISS female connection
64-0604	2800 cc (1 gallon) Polycarbonate bottle and cap assembly with DISS female connection
64-0605	1200 cc (1/2 gallon) bottle and cap assembly mounting basket
64-0606A	Bottle and cap assembly mounting bracket
64-06CAP-6B	Suction cap assembly with suction tubing nipple
64-06CAP	Suction cap assembly with DISS wrench-tight female nut & nipple